VETGuardian PLUS® QUICK START GUIDE

ZOMEDICA®

Set up each session for reliable vital signs.

For full details, see the VETGuardian PLUS User Manual.

START A RUN

Eight steps to begin a monitoring session

1. Power on - Wait for the green Wi-Fi light
2. Log in - Go to myzomedica.com and sign in to your account.
3. Select VETGuardian Plus - Choose the tile for the correct monitor.
4. Choose patient - Search or add; Use the pencil icon to review and edit patient information as appropriate and Save.
5. Center patient - Center the Aim feature brackets on the chest, while keeping the patient's face in the thermal view.
6. Start session - Select Start.
7. Calibrate temperature - For best accuracy, add a manual temperature.
8. Allow the device time (30 - 60 seconds) to begin displaying vitals.

BEST PRACTICES

Positioning and conditions that improve results

- **For best results, the patient should be still. VETGuardian Plus is intended for patients who are stationary.**
- Center the **Aim brackets** on the chest. See the positioning guide on the back for optimal placement.
- **Remove interference.** Keep equipment out of view; the patient should be the hottest object in frame.
- **Avoid heated pads and glass cages** for accurate temperature. Heating devices and glass interfere with temperature readings. Use safely. Keep outside oxygen chambers

NEED HELP?

Customer Support: 734-369-2555 | customersuccess@zomedica.com | Monday-Friday, 8:00 a.m.-5:00 p.m. ET

For installation, troubleshooting, and detailed operating guidance, see the VETGuardian PLUS User Manual.


Visit myZomedica/resources for support tools

© 2026 Zomedica Inc. and Zomedica Inc. All rights reserved. ZOMEDICA and VETGuardian PLUS are trademarks of Zomedica Inc.

Quick reference only - follow the complete Instructions for Use for all operating and safety information.


Proper Positioning

- Position the monitor so it is angled toward the patient's chest. A **45º angle to the animal** is recommended for optimal readings.
- Center the **Aim brackets** on the patient's chest. When practical, keep the patient's face in view.

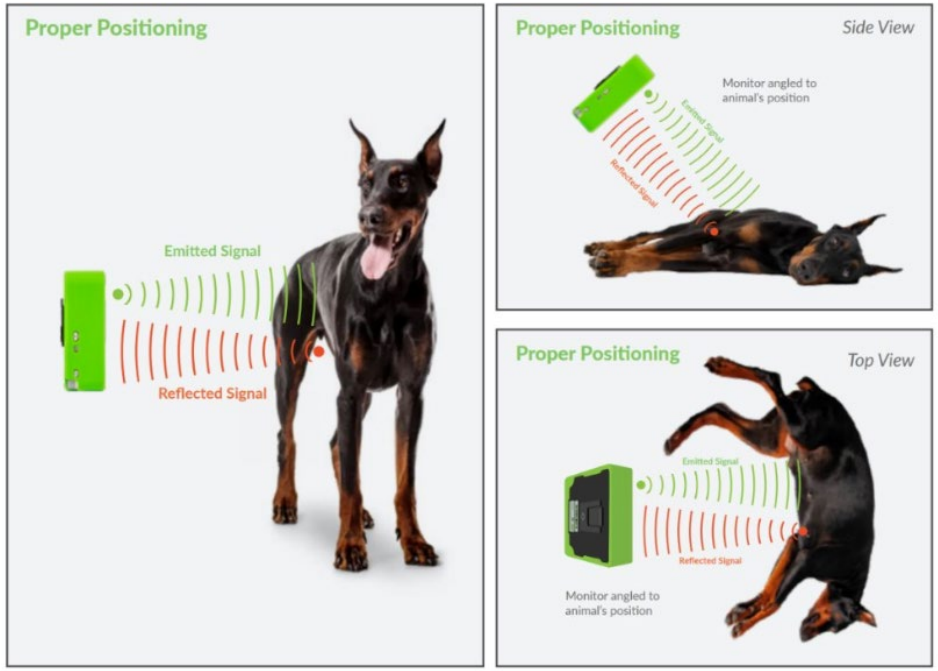



Common Mistakes

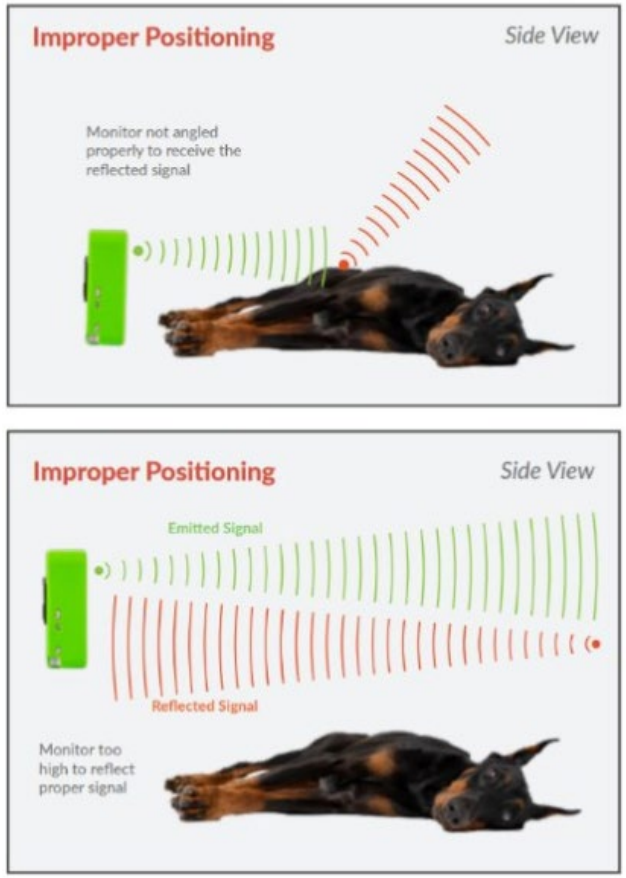

- Reposition the mount if the animal moves and the **Aim brackets** are no longer centered on the patient's chest.
- For accurate temperature readings, avoid heated blankets or other heat-generating items.